




g-Cath™ Tissue Anchor Delivery Catheter with Expandable Tissue Anchors™



g-Cath™ Tissue Anchor Delivery Catheter

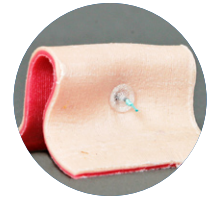

The **g-Cath™ Tissue Anchor Delivery Catheter**

with Expandable Tissue Anchors™ is an integral part of the g-Prox™ tissue approximation system.

The catheter ends in a hollow needle and contains a pair of Expandable Tissue Anchors™. A sliding button on the g-Cath™ allows the surgeon to manually advance the tissue anchors and deliver them on either side of a tissue fold or apposition. Then, the operator uses the g-Cath to bring the two anchors together with adjustable tensioning (similar to a traditional suture).

The tension flattens out the mesh anchors, safely securing the tissue.

Device Features:

- Expandable Tissue Anchors™ spread force for excellent tissue healing
- Constructed of polyester suture and nitinol
- Adjustable tensioning
- High strength
- Creates plications or approximations
- Excellent in-vivo durability
- FDA 510(k) cleared



USGI Medical
1140 Calle Cordillera
San Clemente, CA 92673

Customer Service
1-866-788-USGI (8744)
9:00 AM TO 5:00 PM Pacific Time